



MedITEX instructions for clinical IT administrators:

Standard MedITEX WebPortal Installation

1. Right-click on the ZIP-file with WebPortal installation, choose **"Extract All"**, and save everything in a folder named **"WebPortal"**.
2. Install the **"fmsoft_unigui_extjs_0.93.0.996_setup.exe"**

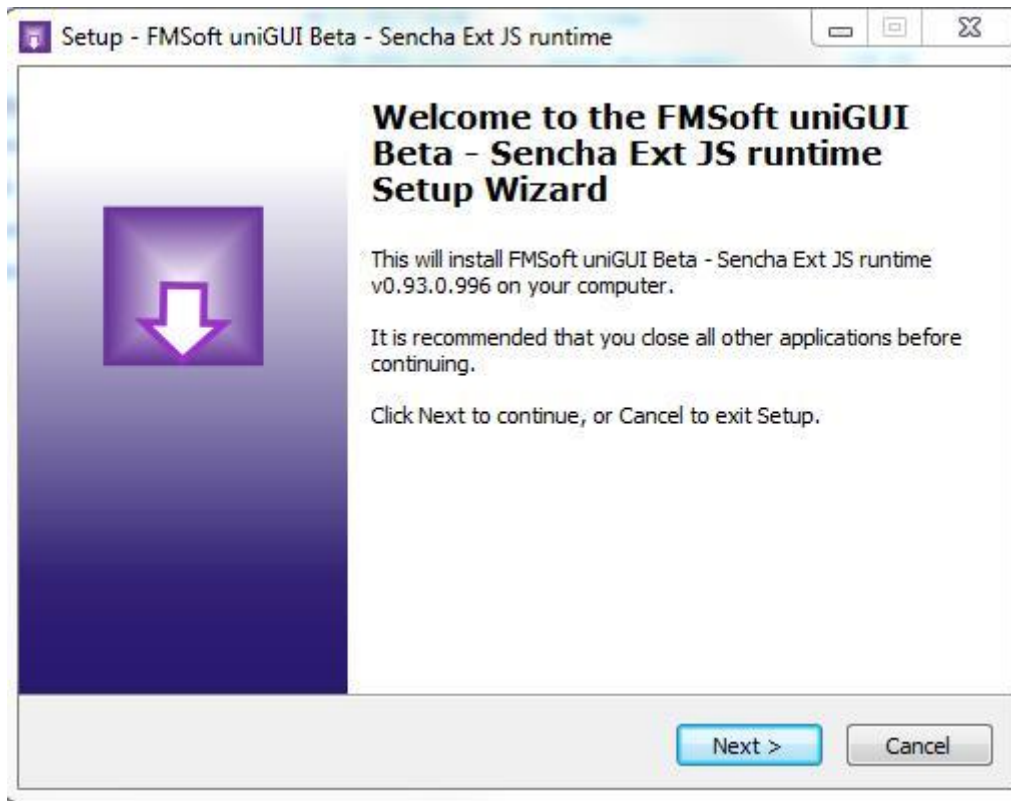


3. Right-click on **"mexWebPortal_Install.bat"** -> **Edit**.
Change this part of the path to the current path of the WebPortal folder.

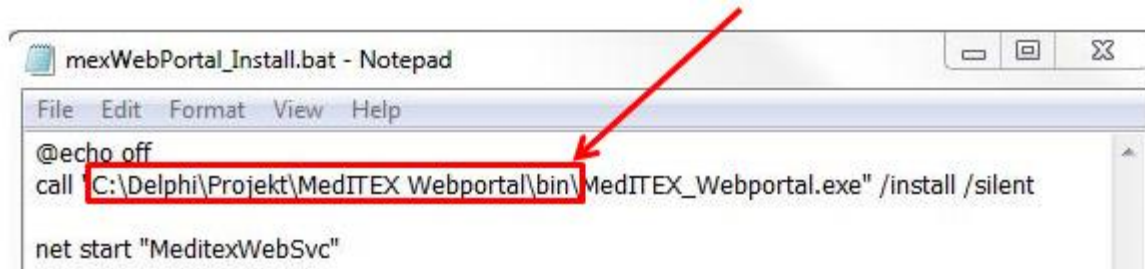

4. Right-click on **"mexWebPortal_Uninstall.bat -> Edit"**.

Change this part of the path to the current path of the WebPortal folder and save it.

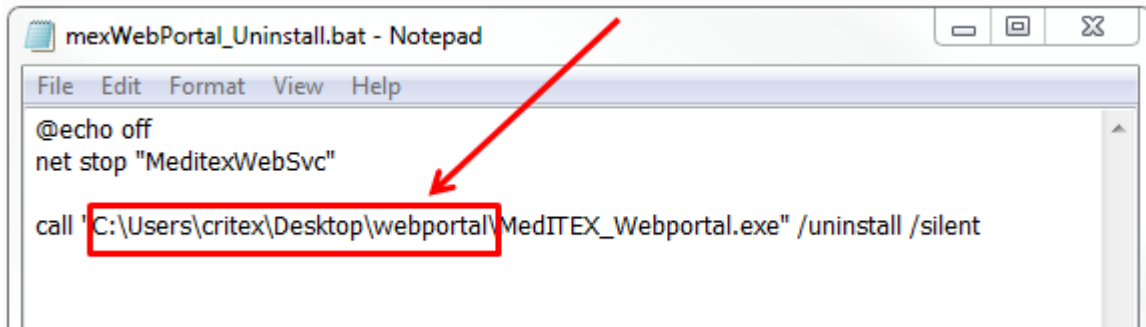


```
@echo off
net stop "MeditexWebSvc"

call 'C:\Users\critex\Desktop\webportal\MedITEX_Webportal.exe' /uninstall /silent
```

5. Right-click on **"Recdate.ini" -> Edit**.

- Change **"localhost"** to the Firebird database location of MedITEX IVF:

There is already a **"Recdate.ini"** file in the **MedITEX VF** folder. Open this file and compare them. If MedITEX WebPortal is not located on the same server as the **IVF**, change **"localhost"** from the WebPortal Recdate.INI file to the same server-id in the **"Recdate.ini"** file from the **IVF** folder. Be sure that the **"DBAlias"** are the same in both **"Recdate.ini"** files.

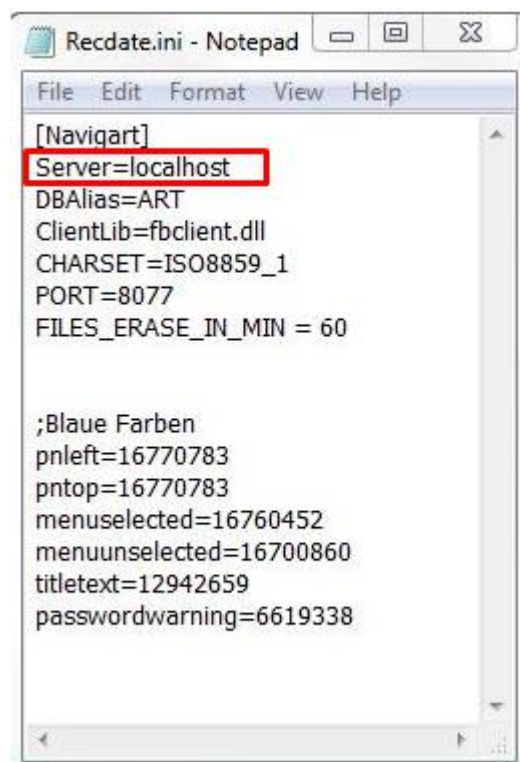


```
[Navigart]
Server=localhost
DBAlias=ART
ClientLib=fbclient.dll
CHARSET=ISO8859_1
PORT=8077
FILES_ERASE_IN_MIN = 60

;Blaue Farben
pnleft=16770783
pntop=16770783
menuselected=16760452
menuunselected=16700860
titletext=12942659
passwordwarning=6619338
```


6. Right-click on **"mexWebPortal_Install.bat"** -> **Run as administrator**

After finishing you should see this window:

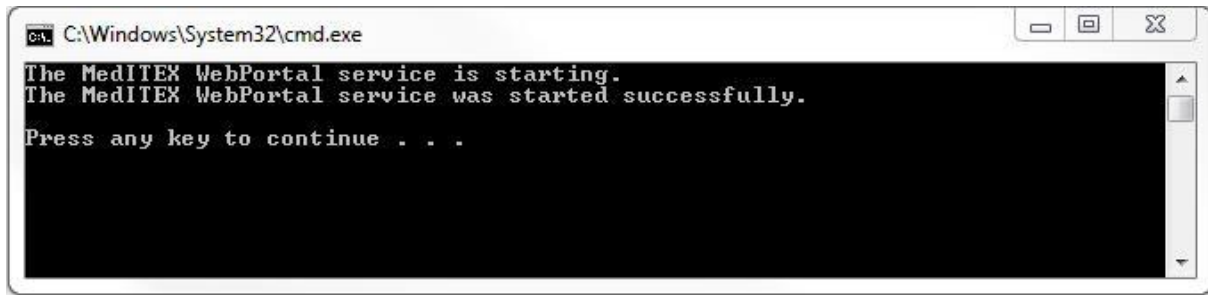


7. Open the **"Control Panel -> Windows Firewall -> Advanced settings"**


Click on **"Inbound Rules"** and choose **"New Rule..."**

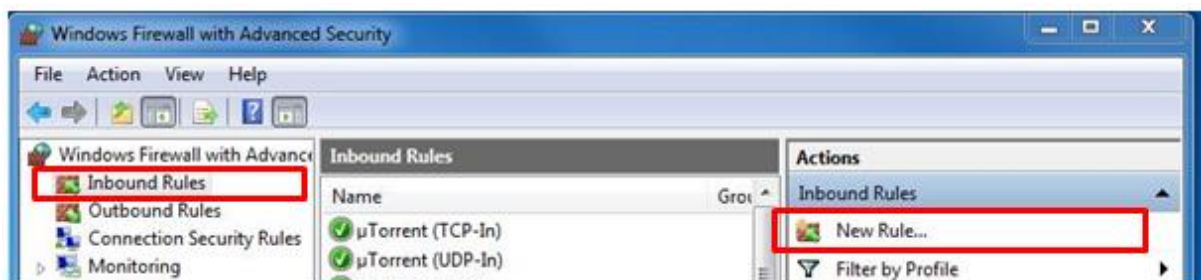

Select **"Port"** and click **"Next"**

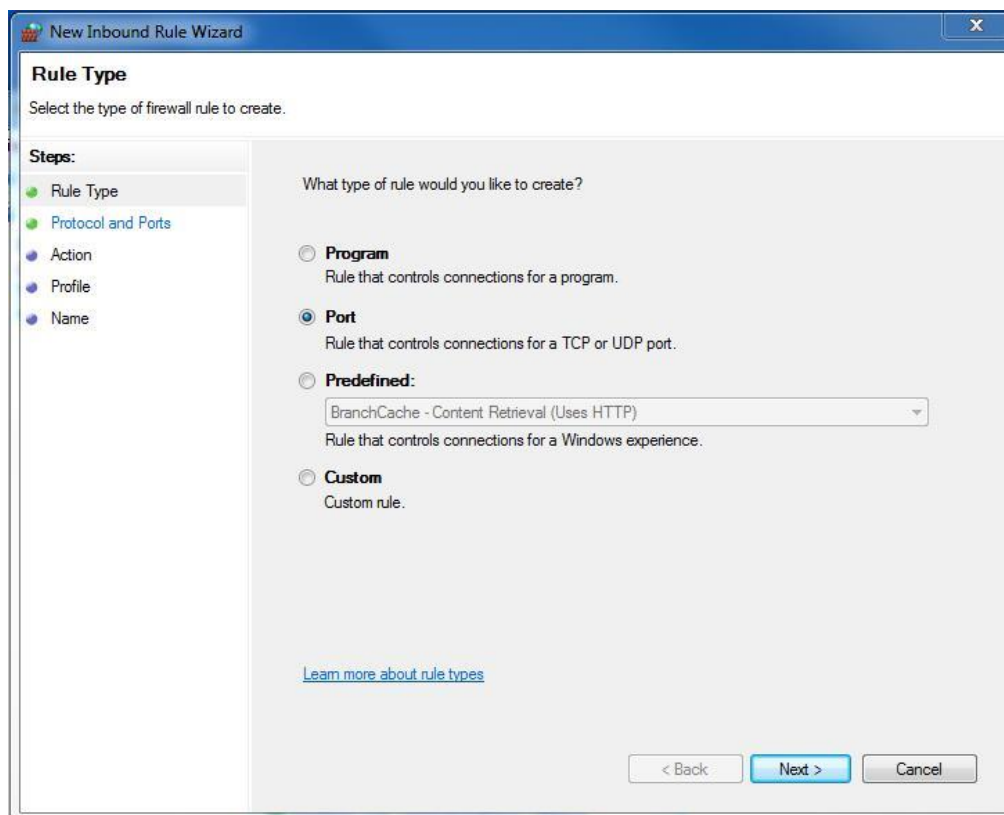


Choose **"TCP"** and **"Specific local ports"** and set the port to **8077** (If this port has already been possessed, you should enter another one – do not forget to change it in RecDate.ini also!)

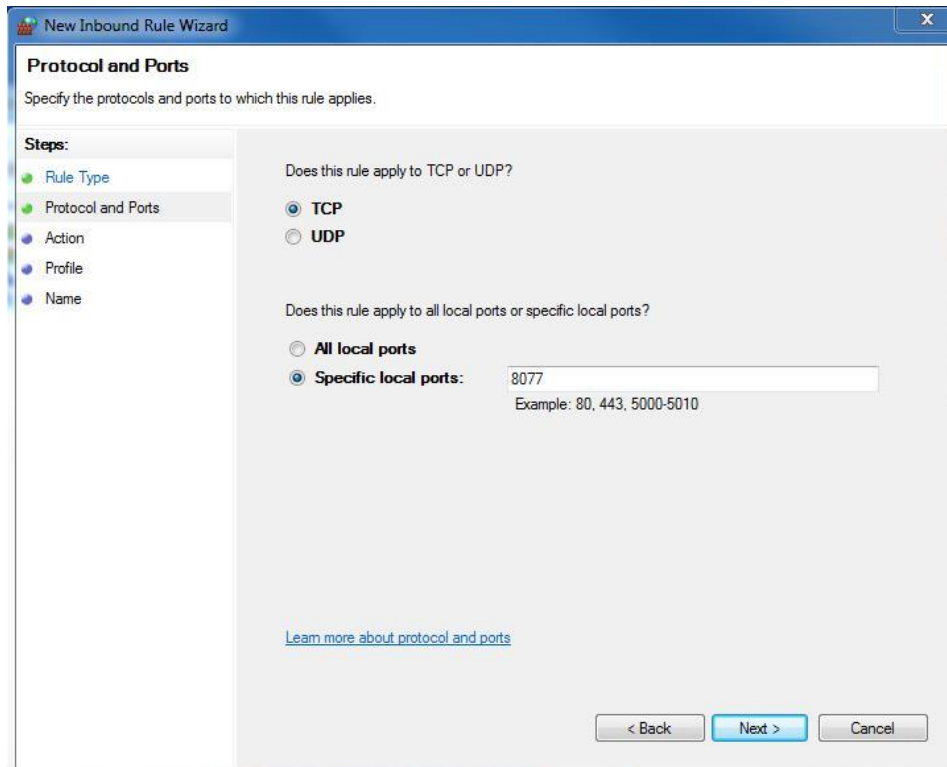

Click **"Next"** until the name is asked for.

Specify the name and optional a description and click **"Finish"**.

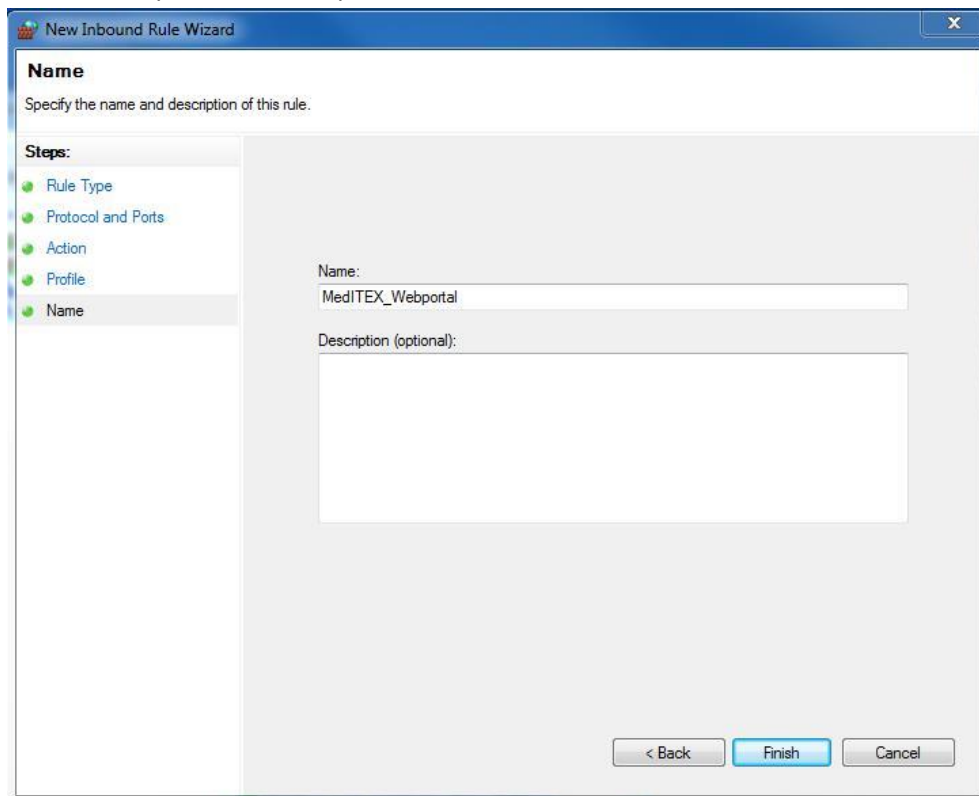

Repeat step 7 with **“Outbound Rules”**:



8. Now it is possible to access the MedITEX WebPortal in a browser with following URL:

http://server-id:8077

Your server-id: The current server-id from step 5 that is also in the **“Reccdate.ini”**.