


## MedITEX instructions for clinical IT administrators: Server Installation Manual for MedITEX IVF

First it is necessary to check if the Firebird ODBC driver and Firebird server are installed on the server. If so, we suggest de-installing the existing Firebird components and replacing them with the ones supplied by us.

### The installation steps for a complete MedITEX IVF server installation:

After you insert the MedITEX IVF setup CD into CD/DVD drive the installation program should start automatically. If the installation does not start, open the contents of the CD and double-click on **“MedITEX IVF.exe”**.

To start the server installation please select **“Server Installation”** option in the installation dialog 4 and press the button **“Next”**.

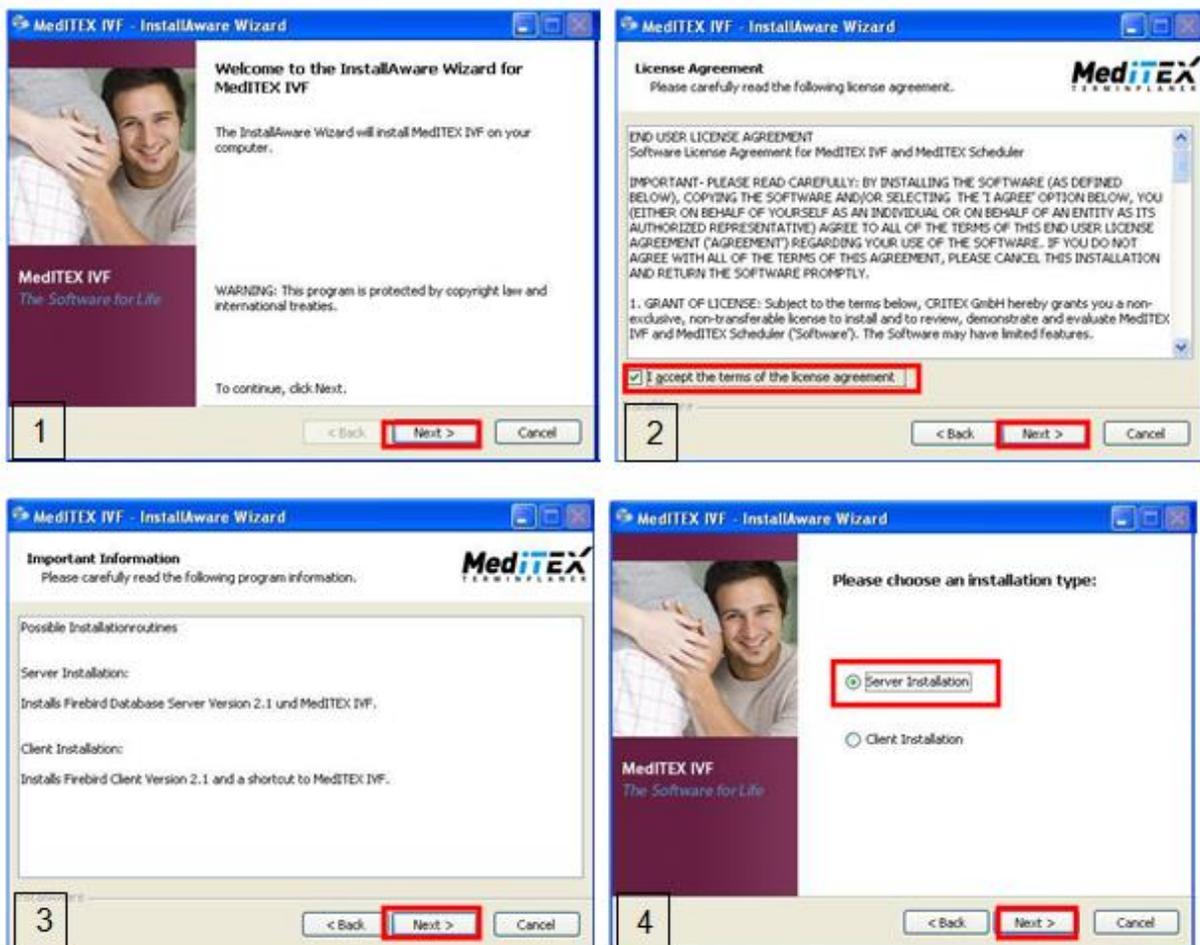

In the following step you have to define the installation path for MedITEX IVF.



**We strongly recommend to install the software in the suggested path "`<drivename>:\Critex`". In this way will the software installation and future maintenance be standardized and simplified.**


**We also recommend not to change the suggested paths during the installation process and accept the offered ones. In this way will the software installation and future maintenance be standardized and simplified.**

### **Installation of Firebird Server Components**


After you have accepted the license terms, please confirm the offered path for the firebird installation.




The configuration of the Classic Server should be set up as in the picture below:

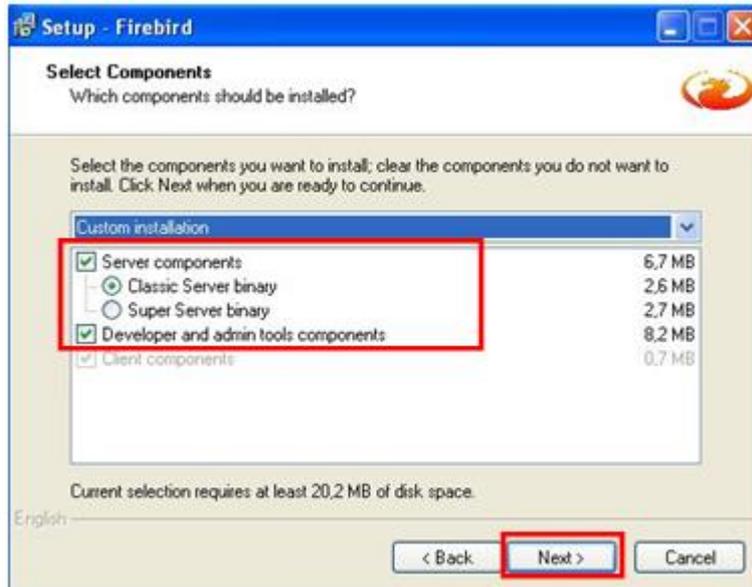


In the next step please select the option „*Don't create a Start Menu folder*“ and click „*Next*“.

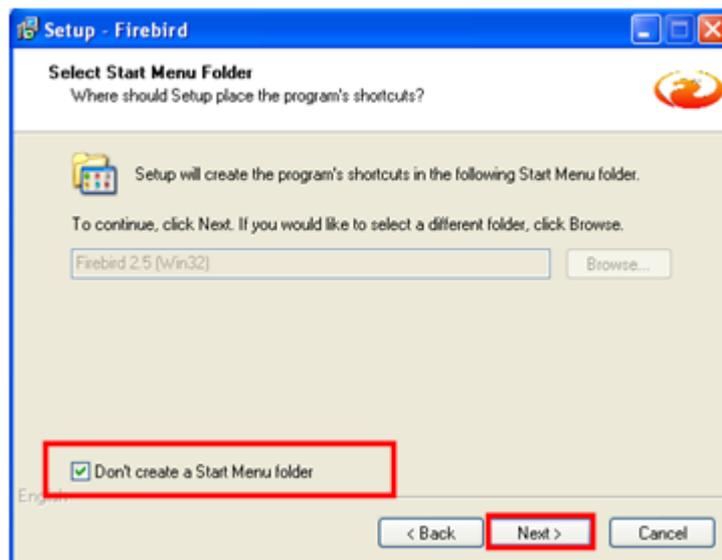

In the next window please set the option „*Copy Firebird client library to <system> directory*“ and click „*Next*“:

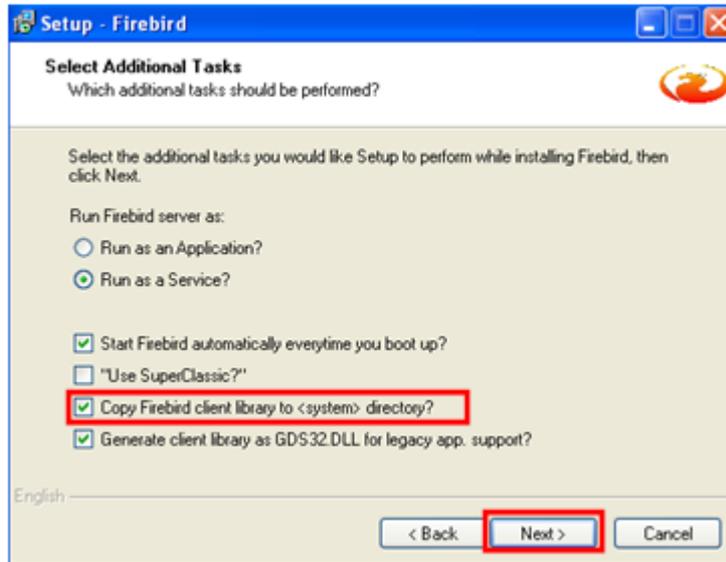


Finally click on „**Install**“:

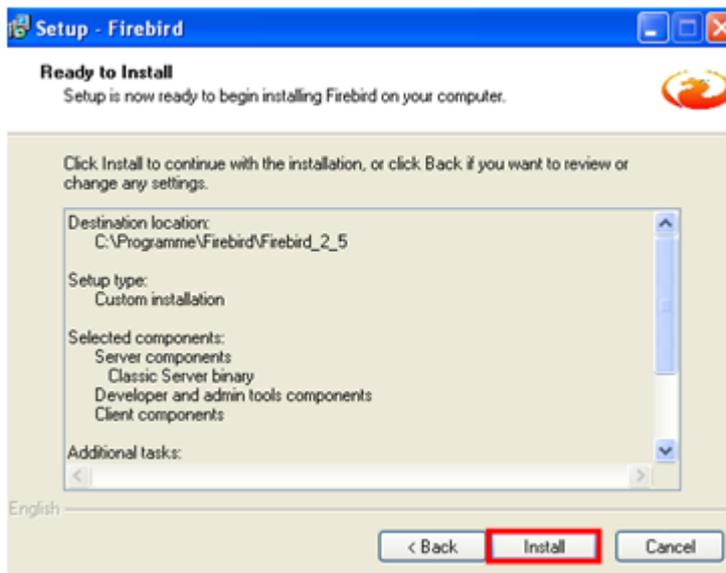

In the last step please empty the check box „**After installation - What Next?** “ and finish the installation by clicking on button „**Finish**“.



Finally the MedITEX IVF Components will be automatically installed.  
After the installation was successfully completed the following message will appear:

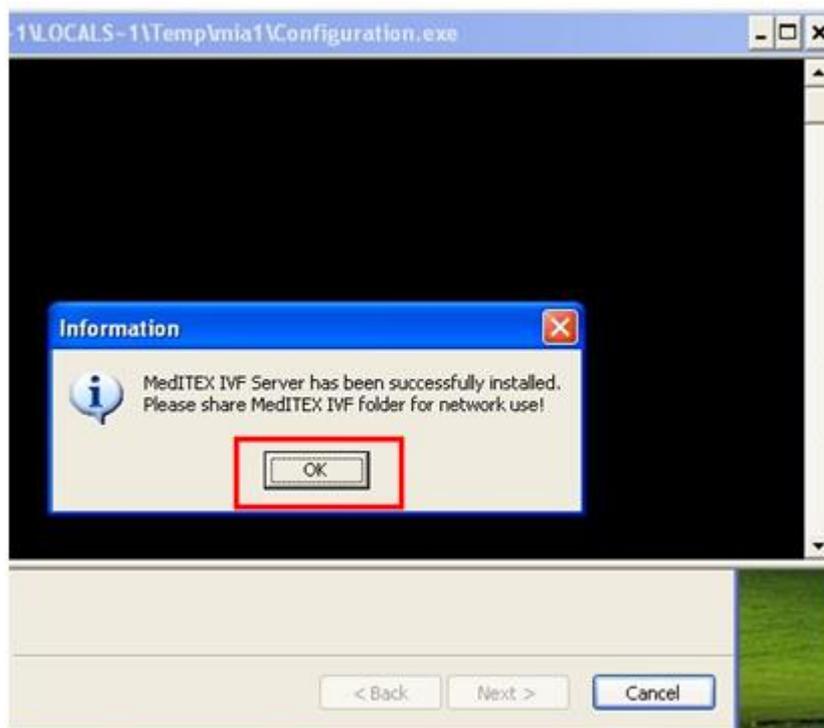


After clicking on the „**OK**“ button the following dialog box will appear.  
Please empty the check box „**Run MedITEX IVF now**“ and click on „**Finish**“.




In the last step it is necessary to share the MedITEX IVF folder for all users. The users need all the rights (read, write and execute) on this folder.

The shared folder is the internal **“MedITEX IVF”** folder. If the software was installed in **“C:\Critex”** it is necessary to share the folder **“C:\Critex\MedITEX IVF”** instead.

In the following pictures we show which folder configurations must be shared, when the previous structure was selected:

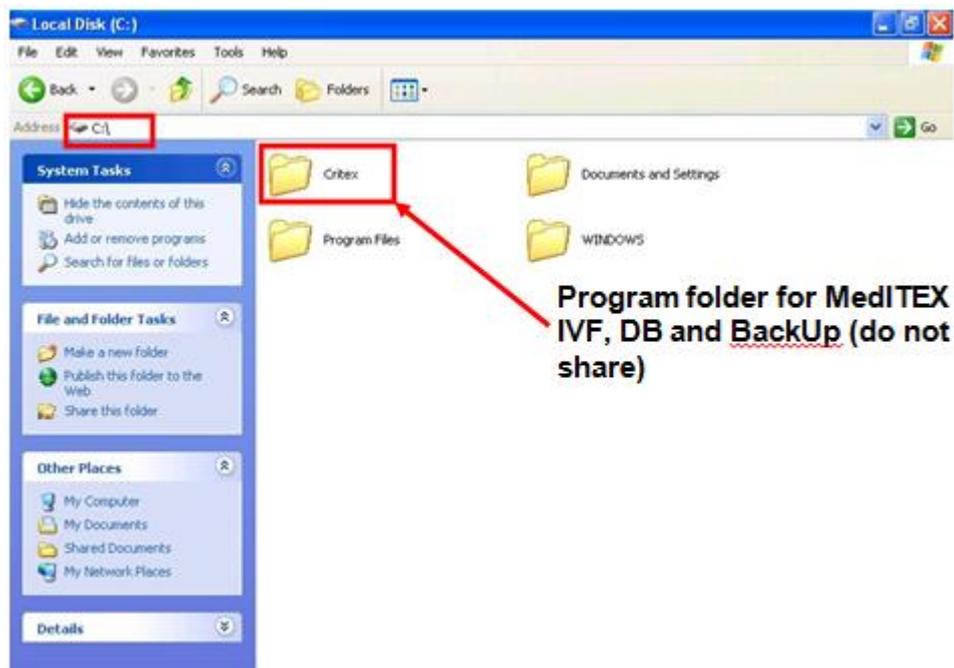

The MedITEX IVF folder must be shared under the name **„MedITEX IVF“**.

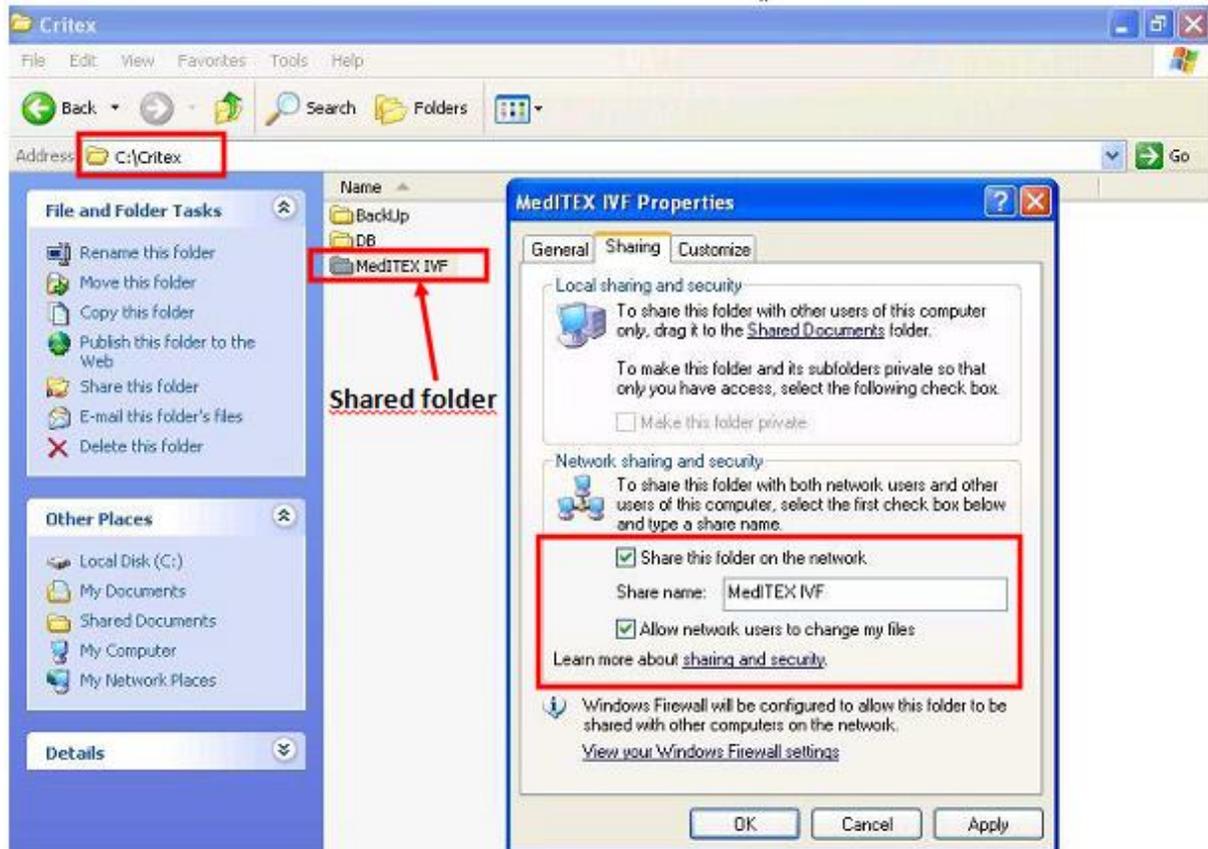

After the successful folder sharing is the installation completed.